

HygroLab
Bench Top Humidity Temperature Indicator
Instruction Manual
V4.0



Contents

Overview.....	3
General Description.....	4
Power Requirements.....	4
Probes and Probe Inputs.....	4
Analog Outputs (optional).....	6
Connection to a PC via RS-232 Serial Port.....	6
Configuration and communication software.....	6
RS485 Networking (HygroLab 2 and HygroLab 3).....	7
Unit System.....	7
Display.....	8
Keypad.....	8
Operation.....	9
Probe Input Selection.....	9
Function Menu.....	9
CALCULATE (HygroLab 2 and HygroLab 3).....	10
MODE (HygroLab 3).....	10
DISPLAY (HygroLab 2 and HygroLab 3).....	13
ADJUST M.PT (full calibration against a reference environment).....	14
ADJUST 1PT (1-point adjustment against a reference environment).....	15
ADJUST REF (1-point adjustment against a reference probe).....	16
PROBE.....	17
SETTINGS.....	17
SYS STATUS.....	18
Errors and Status Messages.....	18
Connectors.....	19
Internal Service Connector (HygroLab 1).....	19
Optional Analog Output Connector (HygroLab 2 and HygroLab 3).....	19
RS232 / RS485 (HygroLab 2 and 3).....	20
RS-485 multi-drop network (HygroLab 2 and HygroLab 3).....	20
Environmental Limits.....	21
Maintenance.....	21
Specifications.....	22
Appendix 1: Practical Advice for Measuring Humidity.....	23
Appendix 2: Maintenance of the ROTRONIC probes.....	23
Cleaning or Replacing the Dust Filter of the Probe.....	23
Periodic Calibration Check of the Probes.....	23
Appendix 3: Calibration Basics.....	24
Temperature Calibration.....	24
Humidity Calibration.....	24
Appendix 4: Humidity Definitions.....	26
Appendix 5: Dew Point Accuracy.....	29
Appendix 6: RS232 Communication Protocol.....	30
Appendix 7: Water Activity Measurement.....	33
Water Activity: definition and applications.....	33
Using the HygroLab 2 or 3 to Measure Water Activity.....	35
General Recommendations.....	37
Using the AwE Mode (HygroLab 3).....	38
Using the AwQuick Mode (HygroLab 3).....	39
Appendix 8: Accessories for the HygroLab.....	42

About the version number: this manual is valid for all instruments with firmware version 4.x, where 4.x can be 4.0, 4.1, etc. (see Function Menu, SYS STATUS). Changes to the last digit of the version number reflect minor firmware changes that do not affect the manner in which the instrument should be operated.

Note: both instrument configuration and adjustment of the HygroClip probes require a PC with the optional ROTRONIC HW3 or HW4 software. Instructions for using the software are not included in this manual. These instructions are shipped separately on the software CD ROM.

Overview

The HygroLab is a bench-top laboratory humidity temperature indicator that can be used with a wide variety of probes to meet specific application requirements. The HygroLab operates with an external 9V AC adapter.

The HygroLab is available in 3 different models:

HygroLab 1: basic indicator

- Accepts up to 2 probes: ROTRONIC HygroClip digital probes or analog probes¹⁾
- Simultaneous display of relative humidity and temperature (one probe at a time)
- Software-based probe calibration (1-point or multi-point)²⁾

- 1) ROTRONIC analog probes only – restrictions apply
- 2) Applies only to the ROTRONIC HygroClip digital probes

HygroLab 2: standard indicator

- Accepts up to 4 probes: ROTRONIC HygroClip digital probes or analog probes¹⁾
- Relative humidity, temperature, dew point, wet bulb, mixing ratio, enthalpy, etc., or user defined calculation such as the difference between temperature and dew point²⁾
- Software-based probe calibration (1-point or multi-point)³⁾
- Possibility of using one of the four probes as a reference to do a 1-point adjustment of the other probes³⁾
- RS232 and RS485 serial ports⁴⁾
- Optional analog outputs (repeat analog probe input signals)

- 1) Third-party pressure probe or ROTRONIC analog probe – restrictions apply
- 2) Uses either a fixed barometric pressure value or the measurements from an analog pressure probe for those parameters that require pressure as a computational input. The fixed pressure value can be changed with the optional HW3 or HW4 software.
- 3) Applies only to the ROTRONIC HygroClip digital probes
- 4) The computed parameter (dew point or other) is not sent to these ports. Use of the HW3 or HW4 software on a PC allows real time computation of the dew point (or other) and facilitates both networking and instrument configuration

HygroLab 3: general purpose indicator with water activity functions

Same as HygroLab 2, except for:

- Accelerated water activity measurement with up to 4 probes (AwQuick mode): allows the measurement of most products in typically 5 minutes
- Full equilibration measurement with up to 4 probes (Standard mode) with automatic detection of equilibrium conditions
- Analog probe inputs with 12-bit resolution

HygroLab 1



HygroLab 2 and 3



General Description

Power Requirements

The HygroLab operates with an external AC adapter providing 9 to 15 VDC, 100 mA. The power receptacle is a 3.5 mm Jack female with DC+ center.

Probes and Probe Inputs

Model	HygroLab 1	HygroLab 2	HygroLab 3
Number of probe inputs	2	4	4
Probe types ¹⁾	ROTRONIC Digital ROTRONIC Analog	ROTRONIC Digital ROTRONIC Analog Third Party Analog	ROTRONIC Digital ROTRONIC Analog Third Party Analog
Number of digital input channels	2	4	4
Number of analog input channels	4	8	8
Analog input A/D resolution	10-bit	10-bit	12-bit

1) Unless otherwise specified when ordering, all probe inputs are factory programmed to accept a HygroClip digital probe. Prior to using any analog probe, the corresponding probe input should be re-programmed. With the HygroLab 2 and 3, this can be done with the HW3 or HW4 software after connecting the RS232 port of the HygroLab to a PC (see separate software manual).

HygroClip Digital Probes

The HygroLab is primarily designed for use with the ROTRONIC HygroClip digital humidity temperature probes. These probes permit to take full advantage of all the features and functions of the HygroLab.


The ROTRONIC HygroClip digital probes are highly accurate and are calibrated entirely by means of software (no adjustment potentiometers). Because calibration and other data are stored in the probe non-volatile memory, the probes are fully interchangeable. When a probe requires calibration or has to be repaired, it can be replaced with another probe in a few seconds.

The ROTRONIC HygroClip digital probes are available in different configurations so as to meet the requirements of each application:



L. 100mm (3.9")
Dia. 15mm (0.6")


HygroClip S
measurement in air
max. 85°C (185°F)

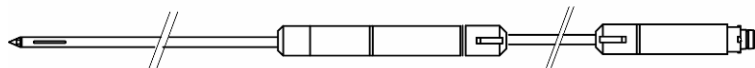

Tube L. 200mm (7.9")
Dia. 5mm (0.2")

Cable: 2m (6.5ft)


HygroClip SP05
measurement in air ducts
max. 85°C (185°F)

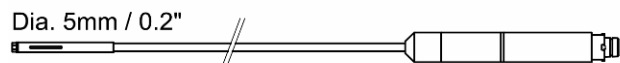

Dia. 5mm / 0.2"

Cable: 2m (6.5 ft)

HygroClip SC05
measurement in tight spaces
max. 100°C (212°F)

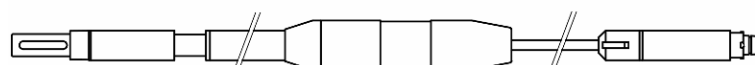

Tube L. 250/400mm/(9.8/15.7")
Dia. 15mm (0.6")

Cable: 2m (6.5ft)


HygroClip HK25 / HK40 air probe
measurement at high temperature
HK25: max. 100°C (212°F)
HK40: max. 200°C (392°F)
wire mesh filter



Tube L. 280mm (11")
Dia. 10mm (0.4")

Cable: 2m (6.5ft)

HygroClip HP28 insertion probe
measurement of materials in bulk
max. 85°C (185°F)

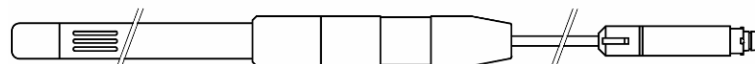


Blade L. 280mm (11")
W. 18mm (0.7")
T. 4mm (0.16")

Cable: 2m (6.5ft)

HygroClip HS28 sword probe
measurement of paper stacks/rolls
max. 85°C (185°F)

Note: all of the above probes are designed for use with the HygroPalm indicator and have a DAT05 connector. An adapter cable MOK-01-B5 is required for connection to the HygroLab (see Appendix 8: Accessories). Similarly, the HygroClip probes designed for use with the HygroFlex transmitter have a T7 connector (not shown here). These probes can be connected to the HygroLab with the T7-03-B5 adapter cable.

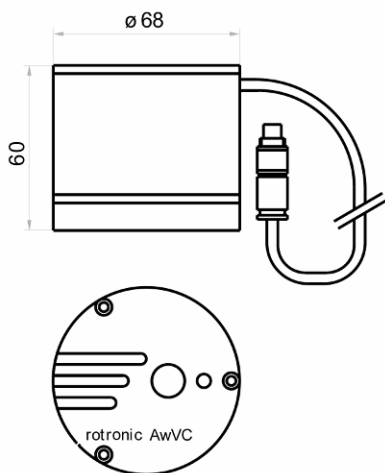

AW-DIO and AwVC-DIO water activity probes

Direct connection to the HygroLab
(probe has a B5 connector)

Analog Probes

The HygroLab can be used with one or more analog probes. The corresponding probe input(s) as well as the scale and unit of the analog input signal(s) should be programmed with the HW3 or HW4 software prior to using the probe(s). For example, an analog pressure probe may be used to provide the local value of barometric pressure for the computation of parameters such as the wet bulb temperature, mixing ratio or enthalpy.

Use of analog probes is subject to the following restrictions and limitations:

- a. ROTRONIC analog humidity-temperature probes with the standard temperature output signal of $-0.5 \dots 2.0 \text{ V} = -50 \dots 200^\circ\text{C}$: because the HygroLab cannot read a negative voltage signal, temperature measurement is generally limited to values above freezing.
- b. Third-party analog probes: single channel probes (one signal), output signal within the range of 0 to 2.5 VDC, supply voltage: 9 VDC or less, maximum current consumption: 10 mA.
- c. Resolution is limited by the 10-bit A/D converter ¹⁾ used for the analog inputs. This converter provides a theoretical maximum of 1024 counts for an input voltage span of 2.5 VDC. In theory, this gives a resolution of $2.5 / 1024 = 0.00244 \text{ V}$. In practice, it is not possible to get 100% of the counts from an A/D converter and the actual resolution should be about 0.0027 V (typical). For example, if a probe with a temperature signal of $0 \dots 1 \text{ V} = 0 \dots 100^\circ\text{C}$ is being used, the signal resolution will be about $100 \times 0.0027 = 0.27^\circ\text{C}$.

1) HygroLab 3: 12-bit A/D converter

Analog Outputs (optional)

Analog outputs are optional on the HygroLab 2 and 3. These outputs simply repeat the analog inputs (if any analog probe is being used). See Connectors – Analog Outputs for details.

Connection to a PC via RS-232 Serial Port

The HygroLab 1 has an internal RS-232 service connector. Both the HygroLab 2 and 3 have also an external RS-232 connector (located on the back panel). Operations such as instrument configuration and probe adjustment require a PC with either the ROTRONIC HW3 or HW4 software.

- HygroLab 1: use of the RS232 port (internal service connector) is limited to configuring the instrument (unit system, output range, etc.). Measurement data and calibration functions are not available on this port.
- HygroLab 2 and 3: use the RS232 port to configure the transmitter, adjust the probe and connect the transmitter to a network. Both instrument configuration and probe adjustment require either the ROTRONIC HW3 or HW4 software. In principle, any communication software can be used to read measurement data from the HygroLab with the commands described in appendix 6 (RS-232 communication protocol).

Configuration and communication software

Both the ROTRONIC HW3 and HW4 software (optional) allow real time computation of the dew point (or other parameter) by the PC while providing additional functionality such as:


- RS-485 networking (HW3: up to 32 units / HW4 or custom software: up to 64 units)
- Instrument configuration (unit system, etc.)
- Adjustment of the ROTRONIC HygroClip digital probes with or without HygroFlex
- Data logging to a PC disk file
- Graphic functions (both on and off-line)

HW3 is compatible with Windows 98 or NT and requires a baud rate of 19200 to be able to communicate with an instrument. HW4 is compatible with Windows 98SE, XP, NT4 with SP6a or higher and allows a baud rate of either 19200 or 57600. For more details see separate instruction manual provided with the software.

RS485 Networking (HygroLab 2 and HygroLab 3)

The optional HW3 or HW4 software offers a simple means of quickly establishing a PC based network.

Up to 32 HygroLab indicators and other devices can be connected together on a network. Any instrument can be used either as a slave or a master, without special configuration. Each unit must be given a unique network address with the HW3 or HW4 software (0 to 31, respectively 0...63). The master is automatically the unit that is connected to the COM port of the PC by means of the RS232 port. See Connectors – RS232 / RS485 for details.

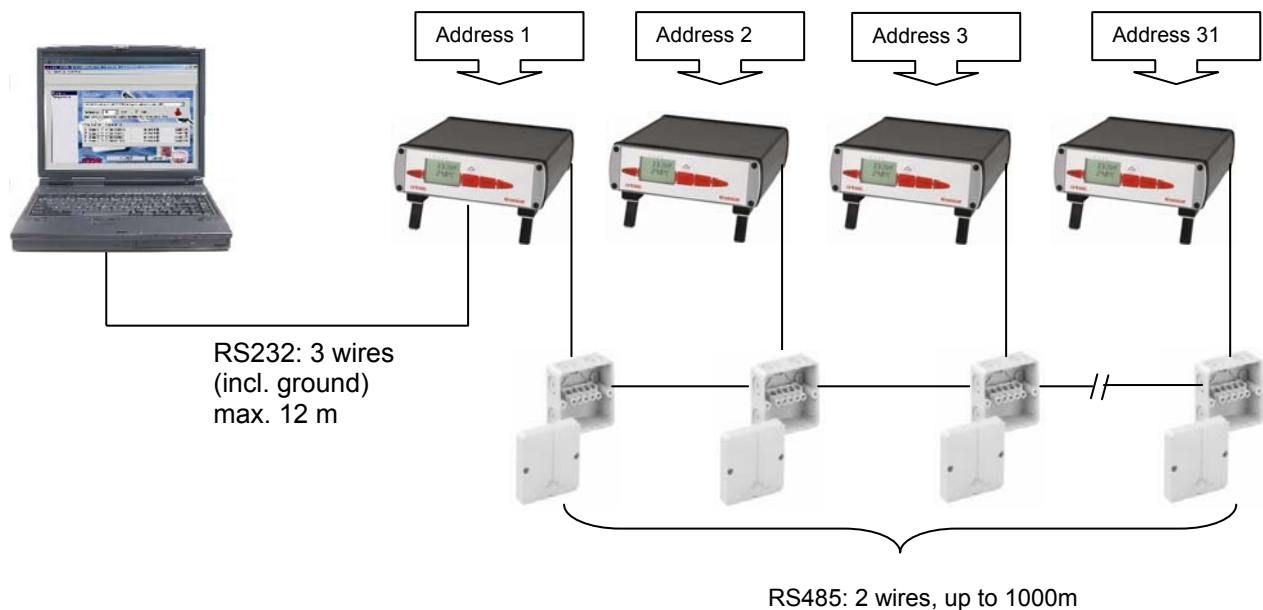

Unit System

The unit system can easily be set to Metric or English from the keypad (except HygroLab 1) -see MENU, SETTINGS. The unit system can also be changed from a PC with the optional HW3 or HW4 software installed. In that case, connect the HygroLab to a free COM port of the PC using either the internal service connector (all models) or the RS232 port (HygroLab 2 and 3)


Note: except for the HygroLab 1, the instrument can be configured to display water activity instead of relative humidity: (instrument configuration with the HW3 or HW4 software):

1.000 Aw = 100.0 %RH


To change the factory configuration for humidity either use the keypad MENU, SETTINGS (except HygroLab 1) or use the optional HW3 or HW4 software installed on a PC and the appropriate connecting cable: cable AC1623 for all models (see Accessories) or cable type RS232 9pol ST-BU for the HygroLab 2 and HygroLab 3 (see Connectors for part number).

Display

The LC display shows which probe input is being displayed (small indicator on top of the display) and up to 2 parameters measured by the probe, with the associated engineering unit. When relevant, the message line provides additional information.

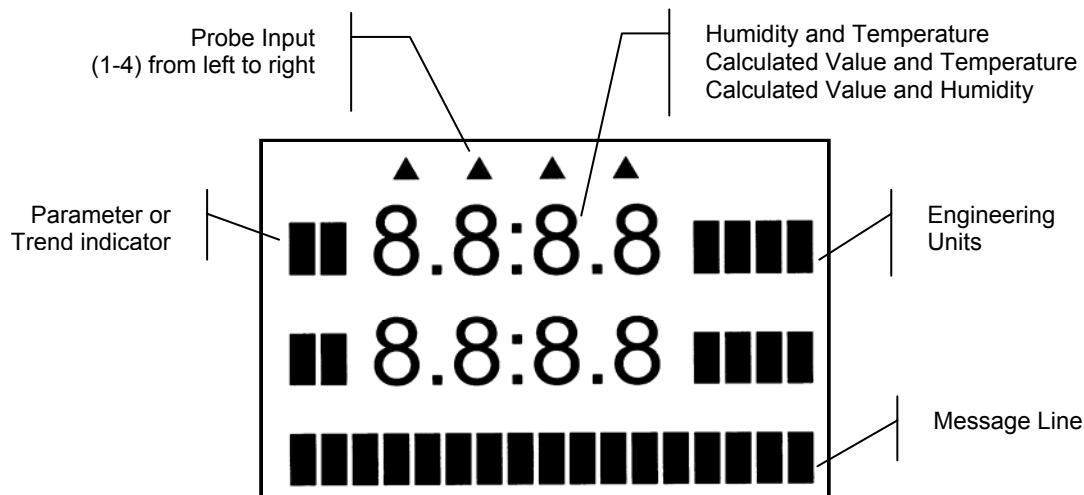

Keypad

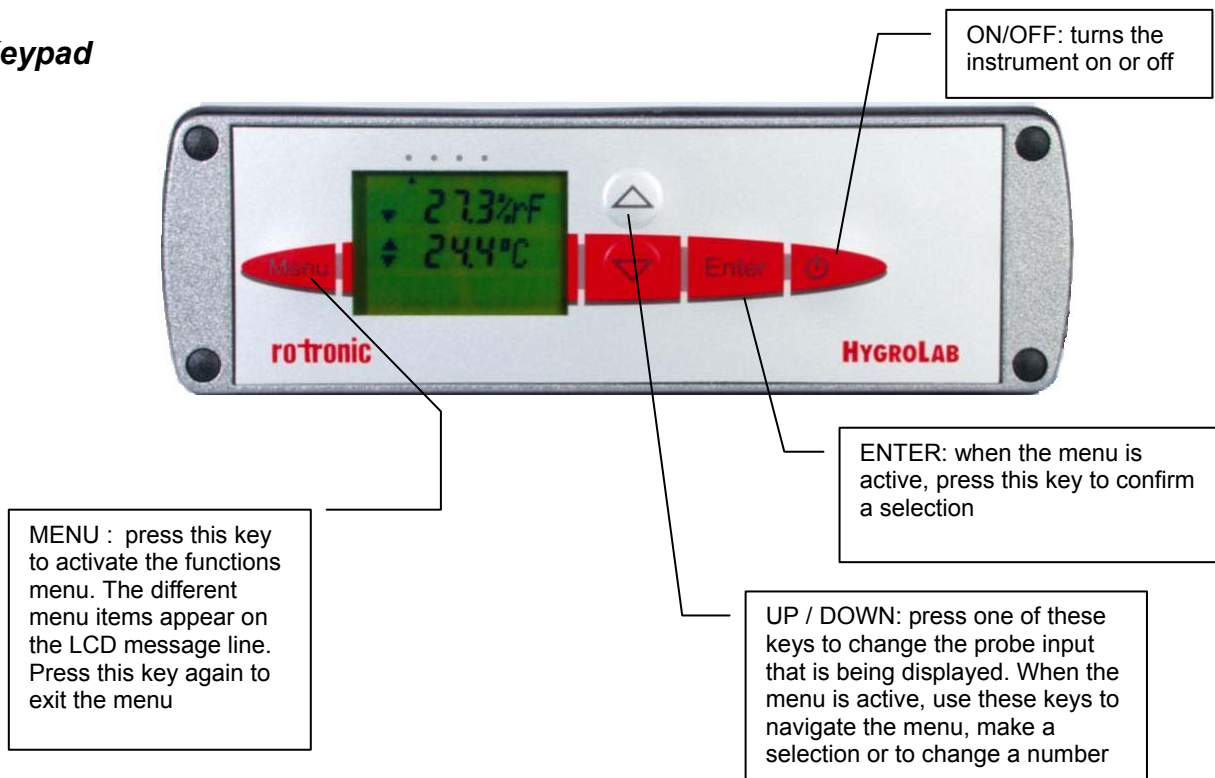

Operation

Connect the AC adapter to the Power connector of the HygroLab and connect a probe ¹⁾. Press on the ON/OFF key. After a brief test and introductory message, the measurements appear on the LC display. Measurements are updated to the display at the rate of 0.75 sec. x number of digital probes (HygroClip) connected to the HygroLab.

1) Unless otherwise specified when ordering, all probe inputs are factory programmed to accept a HygroClip digital probe. Prior to using any analog probe, the corresponding probe input should be re-programmed. This can be done with the HW3 or HW4 software after connecting the RS232 port of the HygroLab to a PC (see separate HW3 or HW4 manual).

HygroLab 1



HygroLab 2 and 3



Probe Input Selection

The display can be switched between probes with the ▼ or ▲ key. The selection is confirmed on the message line of LC display as well as by an indicator located at the top of the display.

Function Menu

To access the function menu, press the MENU key. The first menu item appears on the message line of the LC display. Use the ▼ or ▲ key to navigate the menu. When the desired menu item appears on the message line of the LC display, press the ENTER key to select. Some menu items have sub-items. These can be selected with the ▲ or ▼ and ENTER keys. To exit the menu and return to the normal display mode, press the MENU key. The instrument also returns automatically to the normal display mode when no key is being pressed for some time (main menu: 10 sec., submenu: 30 sec.).

CALCULATE (HygroLab 2 and HygroLab 3)

Definition

This function is used to select or view the humidity parameter that is calculated by the instrument. The calculated parameter is selected individually for each probe connected to the instrument ¹⁾. Prior to entering this function, select on the display the probe to program.

1) Does not apply to third-party analog probes.

The unit system (**Metric** or **English**) can be changed with the keypad (see Functions, Settings) or with the optional HW3 or HW4 software, when the HygroLab is connected to a PC.

Selections

Parameter	M	E	HygroLab 1	HygroLab 2	HygroLab 3
Dew Point / Frost Point ¹⁾	°C	°F	N/A	x	x
Wet Bulb Temperature ²⁾	°C	°F	N/A	x	x
Enthalpy ²⁾	J/g	BTU/lb	N/A	X	X
Vapor Concentration	g/m ³	gr/cuft	N/A	X	X
Specific Humidity ²⁾	g/kg	gr/lb	N/A	X	X
Mixing Ratio ²⁾	g/kg	gr/lb	N/A	X	X
Vapor Concentration at Saturation	g/m ³	gr/cuft	N/A	X	X
Partial Pressure of Water Vapor	hPa	PSI In Hg	N/A	X	X
Saturation Pressure of Water Vapor	hPa	PSI In Hg	N/A	X	X

1) The standard factory setting is frost point for values below freezing. This setting can be changed to dew point with the optional HW3 or HW4 software

2) Calculation of this parameter requires barometric pressure as an input. The fixed value used for barometric pressure can be changed either directly from the keypad - see MENU, SETTINGS - or with the optional HW3 or HW4 software (all models). The HygroLab 2 or 3 can also be configured to accept the input from an analog pressure probe (variable pressure value).

MODE (HygroLab 3)

Definition

This function is used to set the operating mode of the HygroLab 3

Selections

- Standard : normal operating mode (relative humidity, temperature, dew point, etc.)
- Q AwQuick : water activity measurement with the AwQuick algorithm for measurements in typically 5 minutes.
- E AwE : conventional water activity measurement with automatic detection of equilibrium conditions.

Note:

The operating mode selection applies to all probe inputs. Upon start up, the HygroLab 3 remembers in which mode it was before being turned off

● **Standard Mode:**

This mode is the default mode of the HygroLab 3. Upon powering the HygroLab 3, the instrument automatically defaults to the standard mode.

In the standard mode, the HygroLab 3 displays the humidity and temperature measured by the selected probe. Trend indicators appear one minute after powering up the instrument (▲ or ▼ arrow to the left of the water activity and temperature values). These indicators are used to detect equilibrium conditions as required for example during probe calibration. The probe is at equilibrium when both the up and down arrow are displayed at the same time.

● **AwE mode:**

This mode is essentially the same as the standard mode with the difference that the HygroLab 3 automatically detects equilibrium conditions and ends the measurement at that time by freezing the display.

● **AwQuick mode:**

This mode accelerates the measurement of water activity and provides a result in typically 5 minutes. When temperature conditions are stable (both at the product and probe), the measurement obtained with the AwQuick mode is generally within ± 0.005 aw of the measurement that would be obtained by waiting for full equilibration of the product and probe.

Depending on the mode that was selected, proceed as follows

a) **Standard Mode:**

No further action is required after selecting this mode. Simply proceed with the measurements.

b) **AwE Mode:**

Upon selecting this mode, the HygroLab 3 displays the following:

0. 0 0 0 1				
E	H	S	[A w / m i n]

Use the ▼ or ▲ arrow to set the definition of stable humidity. In the above example, the HygroLab will consider humidity to be at equilibrium when the rate of change of the humidity signal is less than 0.0001 Aw per minute. Default value: 0.0001 Aw / min

Press the **ENTER** key to record and confirm this setting. The display now shows the following:

0. 0 1		
E	T	S [° C / m i n]

Use the UP or DOWN arrow to set the definition of stable temperature. In the above example, the HygroLab will consider temperature to be at equilibrium when the rate of change of the temperature signal is less than 0.01 °C per minute. Default value: 0.01 °C / min

Press the **ENTER** key to record and confirm this setting. The AwE mode is now active.

When ready to measure, press the ENTER key to access and start the AwE mode. Press the ENTER key again to stop the AwE mode at any time and return to the Standard mode.

Notes:

- When the HygroLab 3 is in the AwE mode, pressing the ENTER key automatically starts the AwE mode, unless the MENU key was pressed immediately before pressing the ENTER key.
- The ON / OFF key is inoperative in the AwE mode. You must stop the AwE mode in order to be able to turn the power of with the red ON/OFF key.

For instructions on using the AwE mode, see Water Activity Measurement.

c) **AwQuick Mode:**

Upon selecting this mode, the HygroLab 3 displays the following:



The screenshot shows a rectangular display area with a black border. At the top center, the number '4' is displayed. Below it, the text 'Q Dwell [min]' is displayed in a monospaced font.

Use the ▼ or ▲ arrow to set the dwell time. In the above example, the HygroLab 3 will wait 4 minutes before processing the humidity data with the AwQuick algorithm.

Default value: 4 min (recommended value)

Press the **ENTER** key to record and confirm this setting. The display now shows the following:



The screenshot shows a rectangular display area with a black border. At the top center, the number '0.01' is displayed. Below it, the text 'Q TS [°C/min]' is displayed in a monospaced font.

Use the UP or DOWN arrow to set the definition of stable temperature. In the above example, the HygroLab will consider temperature to be stable when the rate of change of the temperature signal is less than 0.01 °C per minute.

Default value: 0.01 °C / min

Press the **ENTER** key to record and confirm this setting. The AwQuick mode is now active.

When ready to measure, press the ENTER key to access and start the AwQuick mode. Press the ENTER key again to stop the AwQuick mode at any time and return to the Standard mode.

Notes:

- When the HygroLab 3 is in the AwQuick mode, pressing the ENTER key automatically starts the AwQuick mode, unless the MENU key was pressed immediately before pressing the ENTER key.
- The ON / OFF key is inoperative in the AwQuick mode. You must stop the AwQuick mode in order to be able to turn the power of with the red ON/OFF key.

For instructions on using the AwQuick mode, see Water Activity Measurement.

DISPLAY (HygroLab 2 and HygroLab 3)

Definition

This function is used to select which parameters the instrument displays. The display mode is specified globally for all probes connected to the instrument.

Selections

- **Measurement**

Relative humidity (in %RH or as water activity: Aw) and temperature or the parameter measured by any third party probe connected to the instrument.

- **Calculated + Temperature**

Calculated parameter (as selected under CALCULATE) and temperature. This selection is not available with third-party probes.

- **Calculated + Humidity**

Calculated parameter (as selected under CALCULATE) and relative humidity. This selection is not available with third-party probes.

- **User Calc**

User defined calculation such as the difference between temperature and dew point of any probe. The display shows only the calculated value. See the separate HW3 or HW4 manual for instructions on how to program the HygroLab 2 or 3 to compute a user specified parameter.

ADJUST M.PT (full calibration against a reference environment)

Definition

The Adjust M.PT function permits the full calibration of any ROTRONIC HygroClip digital probe against a known reference environment. This function is designed to permit calibration at 2 temperature values and at up to 4 relative humidity values.

Selections

The Adjust M.PT function offers the following choices (use the UP and DOWN keys):

RHS: Humidity calibration using the ROTRONIC Humidity Standards ¹⁾
Humidity: Humidity calibration using any suitable reference environment
Temperature: Temperature calibration using any suitable reference environment

¹⁾ For humidity calibration, ROTRONIC offers convenient, certified humidity standards to generate known humidity values (for details, see Appendix 3 - Humidity Calibration).

Procedure

When more than one probe is connected to the instrument, select the probe to be calibrated prior to entering the function (this is the same as the probe that was last displayed).

Regarding the Adjust M.PT function, it is important to observe the following rules:

- a) Always calibrate temperature first (if temperature needs to be calibrated)
- b) When calibrating temperature (2 points), always calibrate at the low value first. The instrument is programmed to use the low temperature value to compute the offset and the high temperature value to compute the gain.

T-low < 40 °C (104°F) : used to compute the calibration offset

T-high ≥ 40 °C (104 °F) : used to compute the calibration gain

For best accuracy, we recommend using a T-low value close to 20°C (68°F).

- c) When calibrating relative humidity (2, 3 or 4 points) with the ROTRONIC humidity standards, always follow the sequence 35 %RH, 80 %RH, 10 %RH or 5 %RH, 0 %RH). When using a reference other than the ROTRONIC humidity standards, use reference conditions that are within the following brackets and observe the sequence:

>25 %RH...≤55 %RH ¹⁾ : used to compute the calibration offset

>55 %RH ¹⁾ : used to compute the calibration gain

>1 %RH...≤25 %RH : sensor linearity adjustment

≤ 1 %RH : sensor linearity adjustment

The HygroLab is programmed to automatically recognize these brackets.

1) For best accuracy, we recommend using values close to 35 %RH and 80 %RH

When the probe is at equilibrium with the reference environment, activate the Adjust M.PT function and make the appropriate selections:

- a) **RHS:** press the ENTER key to select this item. The LC display shows the value read by the probe at the time the ENTER key was pressed. The letters "RHS" are displayed on the message line. Use the ▼ or ▲ key to change the humidity value to match the mean value (at 23°C) written on the certificate supplied with the standard. For a faster change, keep the key pressed down. Press the ENTER key when done.

Note: in the RHS mode, the software automatically compensates for the effect of temperature on the humidity standard. No additional correction is required.

- b) **Humidity or Temperature:** press the ENTER key to select. The LC display shows the value read by the probe at the time the ENTER key was pressed. The word "humidity" or "Temperature" is shown on the message line. Use the ▼ or ▲ key to change the humidity or temperature value to match the value of the reference environment. For a faster change, keep the key pressed down. Press the ENTER key when done.

After pressing the ENTER key, the message "sure?" should appear on the LC display. Press ENTER to confirm. Next, the LC display will confirm that the probe adjustment has been successfully completed. At that time, press ENTER to exit or MENU to return to the function Adjust M.PT and do another calibration point.

Carry on in the same manner for each calibration point, following the sequence described above.

Note: you can exit the function at any time (without calibrating the probe) by pressing the MENU key.

ADJUST 1PT (1-point adjustment against a reference environment)

Definition

The Adjust 1PT function permits to do a 1-point adjustment (temperature or humidity) of any ROTRONIC HygroClip digital probe against a known reference environment. This function is limited to a simple offset adjustment that is applied across the entire measuring range.

Warning: a 1-point adjustment is no substitute for a full calibration (2 or more points). Doing a 1-point adjustment can improve accuracy over a narrow range of conditions and may also be detrimental to accuracy at other conditions.

Selections

The Adjust 1PT function offers the following choices (use the UP and DOWN keys):

RHS: Humidity calibration using the ROTRONIC Humidity Standards ¹⁾
Humidity: Humidity calibration using any suitable reference environment
Temperature: Temperature calibration using any suitable reference environment

1) For humidity calibration, ROTRONIC offers convenient, certified humidity standards that generate known humidity values (for more details, see Appendix 3 - Humidity Calibration).

Procedure

When more than one probe is connected to the instrument, select the probe to be calibrated prior to entering the function (this is the same as the probe that was last displayed).

When the probe is at equilibrium with the reference environment, activate the Adjust 1PT function and make the appropriate selections:

- a) **RHS**: press the ENTER key to select this item. The LC display shows the value read by the probe at the time the ENTER key was pressed. The letters “RHS” are displayed on the message line. Use the ▼ or ▲ key to change the humidity value to match the mean value (at 23°C) written on the certificate supplied with the standard. For a faster change, keep the key pressed down. Press the ENTER key when done.

Note: in the RHS mode, the software automatically compensates for the effect of temperature on the humidity standard. No additional correction is required.

- b) **Humidity or Temperature**¹: press the ENTER key to select. The LC display shows the value read by the probe at the time the ENTER key was pressed. The word “humidity” or “Temperature” is shown on the message line. Use the ▼ or ▲ key to change the humidity or temperature value to the value of the reference environment. For a faster change, keep the key pressed down. Press the ENTER key when done.

¹ 1-point temperature adjustment: limited by the firmware to -20...< 40°C (-4...<104°F)

After pressing the ENTER key, the message “sure?” should appear on the LC display. Press ENTER to confirm. Next, the LC display will confirm that the probe adjustment has been successfully completed. At that time, press ENTER to exit.

Note: you can exit the function at any time (without calibrating the probe) by pressing the MENU key.

ADJUST REF (1-point adjustment against a reference probe)**Definition**

Note: this function is not available with the HygroLab 1

When two or more ROTRONIC HygroClip digital probes are connected to the HygroLab, the Adjust REF function uses probe #1 as a reference to do a 1-point adjustment (both humidity and temperature) of all the other probes simultaneously.

Warning: a 1-point adjustment is no substitute for a full calibration (2 or more points). Doing a 1-point adjustment can improve accuracy over a narrow range of conditions and may also be detrimental to accuracy at other conditions.

Procedure

Expose all probes to the same stable environment. Ventilation of the probes is highly recommended. Observe the readings from all probes for complete equilibration with the reference environment.

Press the ENTER key to enter the function. The display asks you to wait for a short time and then shows the temperature and humidity read by the reference probe at the time the ENTER key was pressed. Note that these values are not updated to the display as long as the function is active. The message line of the display shows REF = Probe 1, indicating that probe #1 is the reference (this cannot be changed).

Press the ENTER key to accept. After pressing the ENTER key, the message “sure?” should appear on the LC display. Press ENTER to confirm. Next, the LC display will confirm that the probe adjustment has been successfully completed. At that time, press ENTER to exit. Note: you can exit the function at any time (without calibrating the probe) by pressing the MENU key.

PROBE

Definition

This function displays the version number and serial number of any ROTRONIC HygroClip digital probe connected to the instrument.

Procedure

Prior to entering the function, select the probe to be displayed (this is the same as the probe that was last displayed).

When entering the function, the version number of the probe is displayed first. Use the ▼ or ▲ key to display the serial number of the probe.

SETTINGS

Definition

This function is used to do the following:

- Turn the trend indicators on or off
- change the symbol used for relative humidity or display humidity as Water Activity (aw)
- adjust the fixed value of barometric pressure used in some psychrometric computations (except HygroLab 1)
- select between the Metric or English unit system

Procedure

Trend: use the ▼ or ▲ key to enable or disable the trend indicators.

The trend of humidity and temperature is shown on the display by an arrow (up or down) to the left of the measured value. Both arrows are shown to indicate stable conditions. Stable conditions are defined as rate of change of less than 0.02 %RH / min or °C / min.

Signal stability is first evaluated after 60 seconds into the measurement and is updated every 30 seconds.

Note: when the trend indicator is enabled, the symbol for the calculated parameter (e.g. Dp for dew point) is no longer displayed to the left of the numerical value.

Humi Unit: use the ▼ or ▲ key to change the symbol used for relative humidity (depends on country and language) or to display humidity as water activity.

Press. Const. (except HygroLab 1) : use the DOWN arrow key to select the digit to be changed. Use the UP arrow key to change the value of the digit.

Units: use the ▼ or ▲ arrow key to toggle between Metric and English.

SYS STATUS

Definition

This function displays the software version of the instrument, the RS232 configuration * (baud rate, parity, bits and stop bit), the network address *, any user defined description for the instrument * and the serial number of the instrument.

This information is for display only. Some parameters such as the network address can be changed with the HW3 or HW4 software, after connecting the instrument to a PC.

* except HygroLab 1

Procedure

When entering the function, the software version number is displayed first. Use the ▼ or ▲ key to display the other data.

Errors and Status Messages


The following is a list of coded messages (101, etc.) that the HygroLab may show on the bottom line of the LC display or via RS232 / RS485.

Errors:		
101	checksum error	the checksum test did not pass during RS-communication.
102	bad command	an unknown command was received
103	disallowed command	a command was received that is reserved for production and service
104	unknown probe input	reference was made to a non-existing analog or digital probe input
105	argument error	error in the number of the arguments in the command or in the value of one of the arguments
106	HygroClip communication error	the HygroClip probe does not answer or is not connected
107	calibration error	the difference between the probe reading and the calibration point is larger than the maximum allowed by the INI command.
108	calibration error (overflow)	internal probe error (or the difference between the probe reading and the calibration value is too large)
110	unknown reference probe	the reference probe is not connected or the reference probe input does not exist
111	Temperature error	During humidity calibration, temperature should be within the limits of 0 and 80°C (32 and 176°F)
Warnings:		
120	no adjustment	calibration
121	No HygroClip probe is connected	
Status:		
130	the probe was adjusted	calibration

Connectors

Internal Service Connector (HygroLab 1)


Open the instrument housing to access the service connector. This 3-pin connector is located on the printed circuit board as shown below:

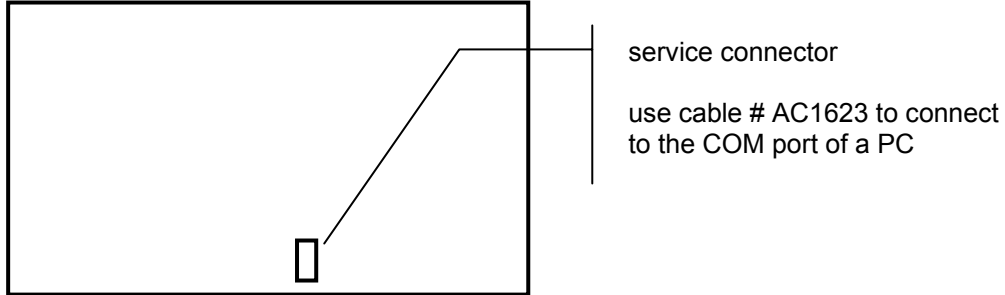

Probe input connectors


Optional Analog Output Connector (HygroLab 2 and HygroLab 3)

If any analog probe is being used, this connector simply repeats the analog probe input signals.

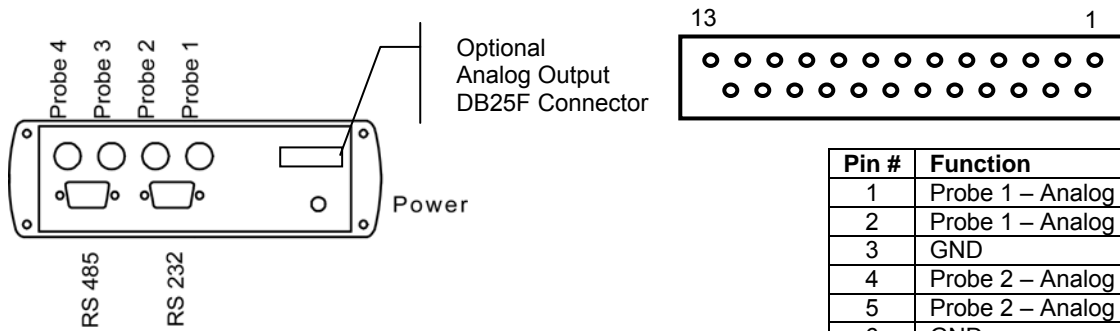


Pin #	Function
1	Probe 1 – Analog 1
2	Probe 1 – Analog 2
3	GND
4	Probe 2 – Analog 1
5	Probe 2 – Analog 2
6	GND
7	Probe 3 – Analog 1
8	Probe 3 – Analog 2
9	GND
10	Probe 4 – Analog 1
11	Probe 4 – Analog 2
12	GND

RS232 / RS485 (HygroLab 2 and 3)

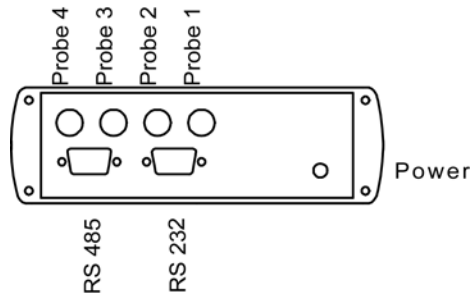


RS232 connector (9-pin Sub D)

Pin #	Function
2	TX
3	RX
5	GND

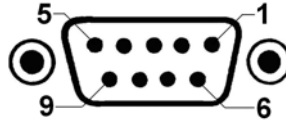

RS485 connector (9-pin Sub D)


Pin #	Function
1	RI +
5	GND (normally not used)
9	RI -

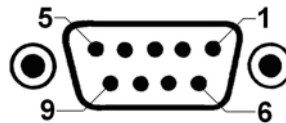

RS-485 multi-drop network (HygroLab 2 and HygroLab 3)

The RS232 connection between the HygroLab and a PC requires 3 wires: RX, TX and GND. Without a signal booster, the maximum cable length is limited to 12 m (39 ft). A connecting cable RS232 9pol ST-BU can be ordered from ROTRONIC (part number 11.01.6218 (1.8 m), 11.01.6230 (3 m), 11.01.6260 (6 m) and 11.01.6290 (10 m)).

The RS485 wiring can be done with a twisted pair (2-wire, telephone type) and wall mounted RJ45 receptacles. No shielded cable and no terminator should be required.

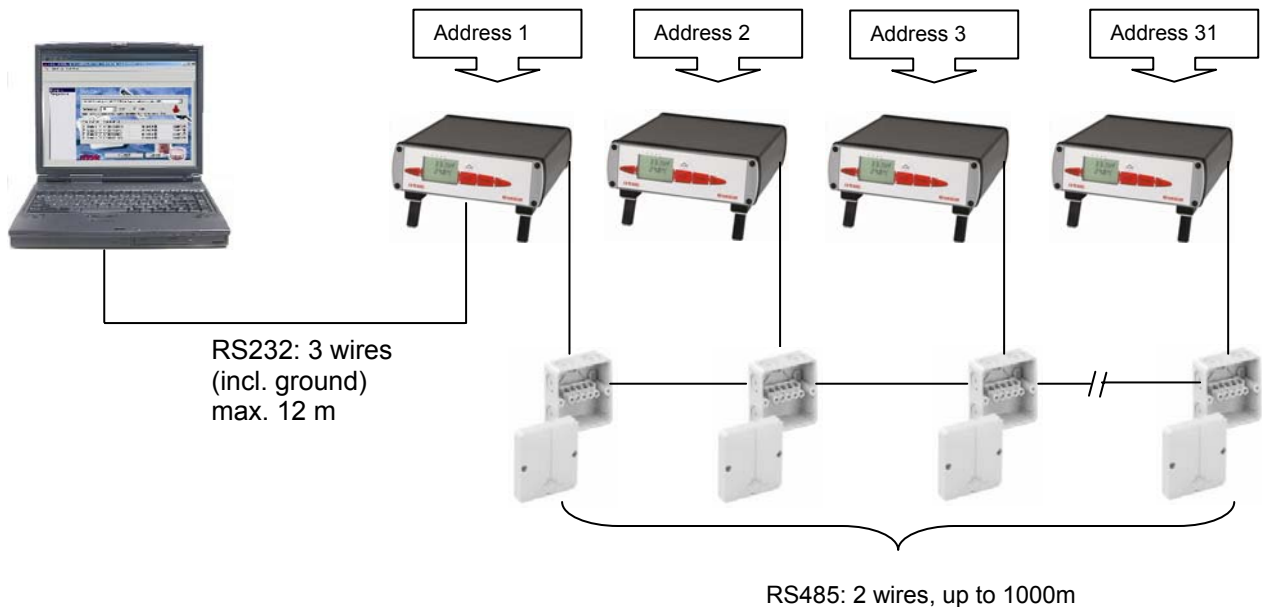

Prior to connecting to the network, configure each individual instrument using either the HW3 or HW4 software (see separate instructions on the software CD ROM). HW3 accepts only a baud rate of 19200. With HW4 a baud rate of either 19200 or 57600 may be used. As far as possible, use a baud rate of 19200 even with HW4. **In any case, be sure to set the baud rate of all devices in the multi-drop to the same value.** Proceed as follows:

- 1) By default, all instruments are shipped by the factory with RS address 0. If you plan on using HW4 for your network, the RS address does not have to be changed as HW4 will do this automatically once the RS-485 network has been established. With HW3, you should connect each instrument, one at a time, to a COM port of the PC and manually configure the RS-485 address so that each instrument has a unique address between 0 and 31. Do the same when you plan on using your own software and use an address range between 0 and 63.
- 2) For easier identification, you may also want to give at that time a unique descriptive name (maximum 30 characters) to each instrument.

Environmental Limits

The HygroLab can operate in the temperature range of $-10...55^{\circ}\text{C}$ ($14...131^{\circ}\text{F}$) and should not be exposed to condensing humidity.

The operating limits at the probe depend on the probe model and are specified separately for each probe.

Maintenance

The HygroLab should not require any routine maintenance. See appendix 2 and appendix 3 regarding the maintenance of the probes.

Specifications

Power supply		AC Adapter 9 to 15VDC, 100 mA DC(+) Tip
Operating limits at electronics		0...99 %RH (non condensing) -10...55°C (14...131°F)
<i>Measured parameters</i> ¹⁾		
relative humidity	indication range	0.0...100.0 %RH
temperature	indication range	-99.9...999.9 °C or °F
barometric pressure ²⁾	indication range	0.000...9999 hPa / PSI / In Hg
<i>Calculated parameters</i> ^{1) 3)}		
dew / frost point ⁴⁾	indication range	-99.9...999.9 °C / °F
wet bulb temperature ⁵⁾	indication range	0.0...100.0 °C / 32.0...212.0 °F
enthalpy ⁵⁾	indication range	-99.9...9999 J/g or BTU/lb
vapor concentration	indication range	0.000...9999 g/m ³ or gr/cuft
specific Humidity ⁵⁾	indication range	0.000...9999 g/kg or gr/lb
mixing ratio ⁵⁾	indication range	0.000...9999 g/kg or gr/lb
vapor concentration at saturation	indication range	0.000...9999 g/m ³ or gr/cuft
part. pressure of water vapor	indication range	0.000...9999 hPa / PSI / In Hg
sat. pressure of water vapor	indication range	0.000...9999 hPa / PSI / In Hg
Number of probe inputs	HygroLab 1	2
	HygroLab 2+3	4
Probe type ⁶⁾	All models	ROTRONIC HygroClip digital probe
	All models	ROTRONIC analog probe ⁷⁾
	HygroLab 2+3	Third-party analog probe ⁷⁾
A/D resolution (analog probes)	HygroLab 1+2	10-bit
A/D resolution (analog probes)	HygroLab 3	12-bit
Display		Alphanumeric LC display
Serial output (except HygroLab 1)		RS 232 / RS485 ⁸⁾
Housing material		Aluminum
Housing dimensions		220 x 170 x 55 mm (8.66 x 6.69 x 2.16")
Weight		1.0 kg (2.2 lb)

1) Specify engineering units when ordering (can be changed by user with optional HW3 or HW4 software). See separate probe specifications regarding the different range limits

2) Not available with HygroLab 1

3) Not available with HygroLab 1. The accuracy of the computed parameter is limited by the accuracy of the measured parameters on which it is based

4) The standard factory setting is frost point for values below freezing. This setting can be changed to dew point with the optional HW3 or HW4 software

5) A fixed pressure value is used for the computations of this parameter. The fixed pressure value can be changed by the user either from the keypad or with the optional HW3 or HW4 software. The HygroLab 2 and 3 can also be configured with the HW3 or HW4 software to accept the input from a pressure probe (variable pressure value)

6) Accuracy, repeatability and operating limits are specified separately for each model of probe

7) ROTRONIC analog probes: because the HygroLab cannot read a negative input voltage, temperature measurement with probes having a standard output is generally limited to values above freezing. Third-party probe (pressure or other): linearized voltage output signal within the range of 0...2.5 VDC, supply voltage 9 VDC, 10 mA max.

8) RS-485 network: requires the first instrument to be connected to the RS-232 COM port of a PC. Compatibility with the ROTRONIC HW3 software requires a baud rate of 19200.

Appendix 1: Practical Advice for Measuring Humidity

The most common source of error when measuring relative humidity is a difference between the temperature of the probe and the temperature of the environment. At a humidity condition of 50 %RH, a temperature difference of 1°C (1.8 °F) typically results in an error of 3 %RH on relative humidity.

When using a humidity probe with an indicator, it is good practice to monitor the display for temperature stability. The probe should be given sufficient time to equilibrate with the environment to be measured. The larger the initial temperature difference between the probe and the environment to be measured, the more time temperature equilibration requires. This time can be shortened, and errors avoided, by using the probe configuration that fits best for your application.

In extreme situations, condensation may occur on the sensors when the probe is colder than the environment. As long as the humidity / temperature limits of the humidity sensor are not exceeded, condensation does not alter the calibration of the sensor. However, the sensor has to dry out before it can provide a valid measurement.


Non-moving air is an excellent insulator. When there is no air movement, surprising differences in temperature and humidity can be noted over short distances. Air movement at the probe generally results in measurements that are both faster and more accurate.

Appendix 2: Maintenance of the ROTRONIC probes

Cleaning or Replacing the Dust Filter of the Probe

Most ROTRONIC probes come with one of the following types of dust filter: (a) protective metal base with a removable filter cartridge or (b) plastic slotted cap with built-in (not removable) filter element.

Depending on the conditions of measurement, the cleanliness of the filter should be checked from time to time. If the probe has a removable filter cartridge, this can be easily removed for cleaning.

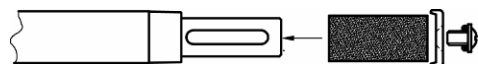

If the probe has a plastic slotted cap with a built-in filter element, cleaning should be done without removing the filter from the probe. In that case, do not use detergents, solvents or other strong chemicals. Either brush the filter or use a little bit of clean water.

Corroded, discolored or clogged filters should be replaced. If the probe has a removable cartridge, simply replace the cartridge (leave the metal base on the probe).

If the probe has a plastic slotted cap with built-in filter element follow these instructions:

- 1) Unscrew the filter from the probe and pull it straight away, in the alignment of the probe, so as not to catch the humidity and temperature sensors.
- 2) Before putting on a new dust filter, check the alignment of both sensors with the probe. The wires that connect the sensors to the probe are very thin and bend easily. If necessary, correct the alignment by tapping the sensor very gently with a smooth object such as a small plastic rod. Do not use sharp pliers or tweezers as this could puncture the sensor and do not pull hard on the sensor.

Periodic Calibration Check of the Probes

Long term stability of the ROTRONIC Hygromer humidity sensor is typically better than 1 %RH per year. For maximum accuracy, calibration of the probe should be verified every 6 to 12 months.

Applications where the probe is exposed to significant pollution may require more frequent verifications.

Both the Pt 100 RTD temperature sensor and associated electronics are very stable and should not require any calibration after the initial factory adjustment.

For routine calibration checks, the probe should be verified at one or two values of humidity. For more details, see Calibration Basics and Functions – ADJUST M.P.T.

For customers who do not want to get involved with field calibrations, our unique probe exchange program offers an attractive alternative. This program takes advantage of the interchangeability of the ROTRONIC HygroClip digital probes. Replacement probes are fully calibrated and are rehabilitated probes with a brand new humidity sensor and filter.

Appendix 3: Calibration Basics

The following choices are available to calibrate the HygroClip probe(s) used with the HygroLab:

- a) Calibration using the HygroLab display and keypad (see Functions).
- b) Calibration using a PC with the optional or HW4 software (see separate instructions for the HW3 or HW4 software) – This choice is available only with the HygroLab 2 and 3, after connecting the HygroLab to the COM port of the PC.
- c) Calibration of the probe alone (removed from the HygroLab), using a PC with the optional HW3 or HW4 software and the MOK-WIN or T7-WIN calibration cable (see separate instructions for the HW3 or HW4 software).

Note: the HygroLab itself should not require any field calibration.

Temperature Calibration

Note: the stability of the Pt100 RTD sensor used to measure temperature is such that temperature calibration in the field is seldom required.

In order to be able to correctly evaluate the accuracy of the temperature measurements provided by the probe, you should be able to meet the following requirements:

- a) Both the probe and a reference thermometer should be ventilated with the same stream of air. Any dust filter used to protect the sensors should be carefully removed from the probe. If the probe has a protective slotted cap, this may be left on the probe.
- b) Air velocity at the sensor should be within the limits of 200 to 500 feet/minute (1 to 2.5 meters/second). Any comparison between two instruments at a velocity under 200 feet/minute may not be valid. Air velocity above 500 feet/minute may damage the unprotected humidity sensor.
- c) The temperature of the air stream should be practically constant.

If you cannot meet the above requirements, you should not attempt to calibrate temperature.

Humidity Calibration

ROTRONIC provides easy-to-use, certified humidity standards for customers who do not have access to a humidity generator. To use these standards, you will need a calibration device that is suitable for your probe.

Calibration Device

The calibration device is a small airtight container that fits on the probe and seals around the humidity sensor. During calibration, a known reference humidity is produced inside the calibration device by means of a humidity standard (usually an aqueous salt solution). To find out which model of calibration device is suitable for your probe, please consult the probe documentation.

Certified Humidity Standards

The ROTRONIC certified standards are available in boxes of 5 glass ampoules of the same value, which can be stored indefinitely. Standards in the range of 5 to 95 %RH are non-saturated aqueous salt solutions that are precisely titrated at our factory for the right concentration. The 0 %RH humidity standard is made of small granules of a highly porous ceramic that have been dried at a high temperature. A Material Safety Data Sheet is available for each standard. Since most standards are a salt solution, parts which have come in contact with the liquid should be cleaned after each use.

Each box of standards comes with a certificate that provides statistical information on the manufacturing batch of the standard. Information on the effect of temperature on each standard is provided on the cover of each box of standard. When calibrating either with the HygroLab or with the HW3 or HW4 software, the effect of temperature on the standards is compensated by the software and no further correction is required. The value of the standards is not affected by altitude.

Instructions for using the Standards

- Install the calibration device on the probe so that the receptacle (or solution holder) is under the probe. Check for a tight fit and remove the receptacle from the calibration device.
- Place one fiber disc (each box of standards includes 5 discs) in the receptacle of the calibration device. The purpose of this disc is to prevent accidental spilling of the solution inside the calibration device or on the humidity sensor.
- Tap the top of the ampoule so that all liquid drops to the bottom of the ampoule. Snap off top and empty contents on fiber disc. Since the ampoule is made of glass, exercise proper caution (gloves, safety glasses) when snapping off the top.
- Put the receptacle back on the calibration device and make sure that the solution does not come in contact with the sensor: The solution inside the calibration device should never be on top of the sensors.
- Allow at least 60 minutes to insure that the calibration device, the solution and the sensor are in a state of equilibrium. This is verified by monitoring the display.
- After adjusting the probe, remove the receptacle from the calibration device. Throw away the wet disc (non reusable). Thoroughly wash and wipe dry the receptacle.

General Recommendations

During calibration, temperature stability is the single most important requirement. If possible, calibrate the probe at room temperature (18 to 25°C). Room temperature should be stable to $\pm 0.25^\circ\text{C}$ or better during the period of time required for each calibration point. Do not calibrate close to an air vent or a heater, in direct exposure to sun rays, etc.

If using a humidity generator to calibrate the probe, make sure that the probe is as fully immersed in the generator as possible to minimize temperature effects.

Instruments Configured to Display Water Activity

Calibrate, the HygroLab exactly as if it were configured to display humidity in %RH. Do not convert the value of the humidity standard to Aw.

Appendix 4: Humidity Definitions

Relative Humidity and Water Activity

Relative humidity is the ratio of two pressures: %RH = $100 \times p/ps$ where p is the actual partial pressure of the water vapor present in the ambient and ps the saturation pressure of water at the temperature of the ambient.

Relative humidity sensors are usually calibrated at normal room temperature (above freezing). Consequently, it is generally accepted that this type of sensor indicates relative humidity with respect to water at all temperatures (including below freezing).

Ice produces a lower vapor pressure than liquid water. Therefore, when ice is present, saturation occurs at a relative humidity of less than 100 %. For instance, a humidity reading of 75 %RH at a temperature of -30°C corresponds to saturation above ice.

Water activity is the same as relative humidity except for the fact that it is expressed in p.u. instead of %RH: $100.0\% \text{RH} = 1.000 \text{ Aw}$.

Dew Point / Frost Point Temperature

The dew point temperature of moist air at the temperature T , pressure P_b and mixing ratio r is the temperature to which air must be cooled in order to be saturated with respect to water (liquid).

The frost point temperature of moist air at temperature T , pressure P_b and mixing ratio r is the temperature to which air must be cooled in order to be saturated with respect to ice.

Wet Bulb Temperature

The wet bulb temperature of moist air at pressure P_b , temperature T and mixing ratio r is the temperature which the air assumes when water is introduced gradually by infinitesimal amounts at the current temperature and evaporated into the air by an adiabatic process at constant pressure until saturation is reached.

Vapor Concentration

The vapor concentration (density of water vapor in a mixture) - or absolute humidity - is defined as the ratio of the mass of water vapor M_v to the volume V occupied by the mixture.

$$D_v = M_v / V, \text{ expressed in grams/m}^3 \text{ or in grains/cu ft}$$

This can be derived as follows from the equation $PV = nRT$:

a) $M_v = n \times m_w$, where :

n = number of moles of water vapor present in the volume V
 m_w = molecular mass of water

b) $D_v = M_v / V = n \times m_w / V = m_w \times p / RT$, where:

$$m_w = 18.016 \text{ gram}$$

$$p = \text{partial pressure of water vapor [Pa]}$$

$$R = 8.31436 \text{ Pa} \times \text{m}^3 / \text{°K} \times \text{mole}$$

$$T = \text{temperature of the gas mixture in } \text{°K}$$

$$D_v [\text{g} / \text{m}^3] = p / 0.4615 \times T$$

$$1 \text{ gr (grain)} = 0.0648 \text{ g (gram)}$$

$$1 \text{ cu ft} = 0.0283168 \text{ m}^3$$

$$D_v [\text{gr} / \text{cu ft}] = 0.437 \times D_v [\text{g} / \text{m}^3]$$

Specific Humidity

The specific humidity (also known as mass concentration or moisture content of moist air) is the ratio of the mass M_v of water vapor to the mass ($M_v + M_a$) of moist air in which the mass of water vapor M_v is contained.

$$Q = M_v / (M_v + M_a)$$

$$Q = p m_w / (p m_w + (P_b - p) m_a)$$

$$m_w = 18.016 \text{ gram}$$

$$m_a = 28.966 \text{ gram}$$

$$p = \text{partial pressure of water vapor [Pa]}$$

$$p_a = \text{partial pressure of dry air [Pa]}$$

$$P_b = \text{total or barometric pressure [Pa]}$$

$$Q [\text{g} / \text{kg}] = 1000 p / (1.6078 P_b - 0.6078 p)$$

$$1 \text{ gr (grain)} = 0.0648 \text{ g (gram)}$$

$$1 \text{ lb} = 0.4535923 \text{ kg}$$

$$Q [\text{gr} / \text{lb}] = 7 \times Q [\text{g} / \text{kg}]$$

Mixing Ratio by Weight

The mixing ratio r of moist air is the ratio of the mass M_v of water vapor to the mass M_a of dry air with which the water vapor is associated:

$$r = M_v / M_a$$

$$M_v = n \times m_w = m_w \times p \times V / RT$$

$$M_a = n \times m_a = m_a \times p_a \times V / RT = m_a \times (P_b - p) / RT, \text{ where:}$$

$$m_w = 18.016 \text{ gram}$$

$$m_a = 28.966 \text{ gram}$$

$$p = \text{partial pressure of water vapor [Pa]}$$

$$p_a = \text{partial pressure of dry air [Pa]}$$

$$P_b = \text{total or barometric pressure [Pa]}$$

$$R = 8.31436 \text{ Pa} \times \text{m}^3 / \text{°K} \times \text{mole}$$

$$T = \text{temperature of the gas mixture in °K}$$

$$r = m_w p / m_a (P_b - p)$$

$$r = 621.97 \times p / (P_b - p) \text{ [g / kg]}$$

$$1 \text{ gr (grain)} = 0.0648 \text{ g (gram)}$$

$$1 \text{ lb} = 0.4535923 \text{ kg}$$

$$r \text{ [gr / lb]} = 7 \times r \text{ [g / kg]}$$

Enthalpy

The enthalpy (or energy content) of moist air at pressure P_b , temperature t (°C) and mixing ratio r (g/kg) is defined by:

$$h \text{ [kJ / kg moist]} = 1.00464 t + 0.001846 r \times t + 2.5 r$$

Note: by convention, the enthalpy of dry air ($r = 0$) at 0°C is equal to zero. Negative values of enthalpy are possible and indicate that the energy content of the air / vapor mixture is less than the energy content of dry air at 0°C

$$1 \text{ lb} = 0.4536 \text{ kg}$$


$$1 \text{ BTU} = 1.05507 \text{ kJ}$$

$$h \text{ [BTU / lb]} = 0.4299 \times h \text{ [kJ / kg]} + 7.68$$

The value 7.68 is added to make the temperature of 0°F the reference for enthalpy expressed in BTU / lb

Appendix 5: Dew Point Accuracy

The HygroLab 2 and 3 use the relative humidity and temperature measurements to compute another parameter such as dew point, mixing ratio, enthalpy, etc. The accuracy of this conversion varies, depending on the humidity and temperature conditions. Typical accuracy for the dew point conversion is provided in the graph below:

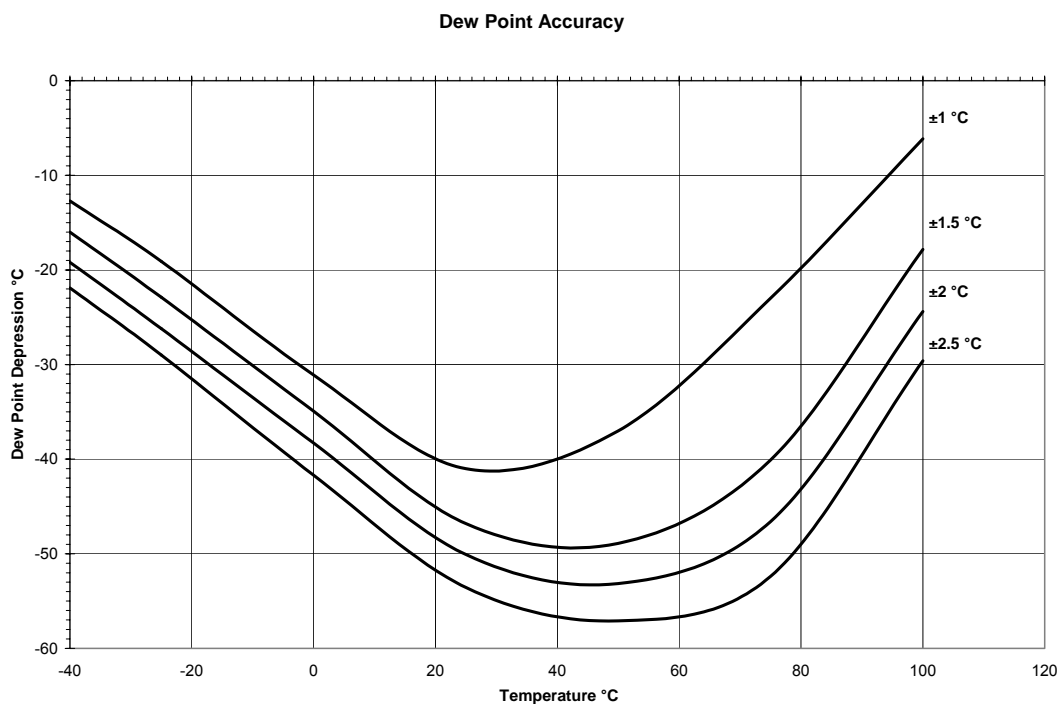

The accuracy of the dew point conversion is less than the accuracy that is normally achieved with a chilled mirror instrument. This is especially true at low dew point values (dew point depression of more than -40 ... -50°C) and at low temperatures. In many applications, repeatability is more critical than accuracy. Repeatability of the dew point conversion is typically 1/3 of the accuracy shown above. Using a dew point conversion permits to measure conditions that cannot be measured with a chilled mirror instrument. Typically, the measuring head of a chilled mirror instrument is limited to a maximum of 70°C and, therefore, cannot measure a dew point temperature above 70°C . In addition to the ability to operate at high temperature, the relative humidity probe offers significant advantages for industrial applications: less maintenance and higher tolerance to contaminants, no sampling system, no uncertainty between dew and frost and better response to fast changing conditions.

Appendix 6: RS232 Communication Protocol

The following describes the RS-232 communication protocol that is used by the instruments from ROTRONIC when connected to the COM port of a PC. This information is provided for users who for any reason do not wish to use the HW3 or HW4 software and for users who plan on using a terminal server (TCP/IP).

Connections

The exchange of data requires 3 wires: TXD (transmit), RXD (receive) and SG (signal ground). When the instrument requires power from the COM port (example: MOK-WIN probe calibration cable), the RTS and DTR wires are also used.

Serial Interface Settings (RS-232 / RS-485)

The HygroLab indicator is designed to accept any of the following serial interface settings:

Settings 1

These settings are the factory default and are used by the ROTRONIC HW3 software. Only 7-bit ASCII characters (ASCII 32 to 127) are being used.

Parity : even
Data bits : 7
Stop bits : 1

When using third-party communications software or a device server (connection to an Ethernet LAN), you should use this group of settings and also set flow control to none. If necessary, force incoming data to 7-bit ASCII.

Settings 2

These settings are used by the ROTRONIC HW4 software.

Parity : none
Data bits : 8
Stop bits : 1

Baud rate:

The HygroLab can communicate with a baud rate of either 19200 (factory default) or 57600.

The ROTRONIC HW4 software (version 1.2.1 or higher) is required to configure the HygroLab for a baud rate of 57600. HW4 is compatible with both 19200 and 57600 baud rates. However, please note that all instruments connected to an RS-485 multi-drop should use the same baud rate.

IMPORTANT: compatibility with the ROTRONIC HW3 software requires a baud rate of 19200. With any other software, we generally recommend using this baud rate because communication with the HygroLab is more reliable than with a baud rate of 57600.

Data Integrity

- 7-bit characters are sent with a parity bit. 8-bit characters are sent without parity bit.
- Strings (made of 7-bit characters) can be sent with or without checksum at the end of the string. When the checksum is not known – for example: when sending a command - the checksum is replaced by the } character at the end of the string.
- For the transmission of 8-bit characters (HW4 software), we make use of predefined initialization and end-of-communication sets of characters with a checksum.

Communication (Serial Interface Settings 1)

Communication with the instrument is always initiated from the PC. Any string sent from the PC should begin with the { character to start the synchronization. This should be followed by a 1-character product identification (see table below). The address is next (this is a 2-digit number between 00 and 31 used to distinguish between the different instruments, probes, functions etc. that are connected to the PC). The address is followed by a command and by any values or parameters that may be required to carry out the command. Use a carriage return (CR) to terminate the command string.

Example: {M03RDD} + CR (asks instrument M03 to send measurement data)

When several instruments are connected together (RS485 network), commands directed to instruments other than the instrument that is physically connected to the COM port of the PC should be preceded with the | character (ASCII 124)

Example: |{U01RDD} + CR

The instrument answers with a string beginning with the { character, the product identification and the address, followed by the command that was received and any data field that may have been requested by the command. This is followed by a checksum.

The length of the string may vary, depending on how much data is included. Since the string itself does not include any information about its own length, a carriage return (CR) is always used to terminate the string.

Example: {M03RDD 0027.50;0023.76;0029.31;0022.77;#M (terminated with CR)

Note: Some instruments may insert the \$ character at the second or third position of the string.

1-character Product Identification (Serial Interface Settings 1)

ID	Product	Data request command
x	MOK-Win	{x00M} + CR
0	HygroPalm 1	Not Applicable
u	HygroPalm 2	{u00RDD}+CR – assuming the address is 0
U	HygroPalm 3	{U00RDD}+CR – assuming the address is 0
1	HygroFlex 1	Not Applicable
m	HygroFlex 2	{m00RDD}+CR – assuming the address is 0
M	HygroFlex 3	{M00RDD}+CR – assuming the address is 0
2	HygroLab 1	Not Applicable
b	HygroLab 2	{b00RDD}+CR – assuming the address is 0

B	HygroLab 3	{B00RDD}+CR – assuming the address is 0
R	M33 series	{R00RDD}+CR – assuming the address is 0
X	Reserved	{X00RDD}+CR – assuming the address is 0

CR=Chr(13)

Special Situations (single instrument – Serial Interface Settings 1)

Using address 99 in a data request causes instruments with any address to answer. For this reason, address 99 should not be used when several instruments are connected to a network. Address 99 can be useful when trying to communicate with a single instrument of unknown address. The address is returned as part of the answer string.

When the product ID is unknown, it can be replaced with a blank in the command. The product ID is returned as part of the answer string.

Example: { 99RDD}+CR

Examples

Data can be read from the instrument directly connected to the COM port of the PC (RS232) as well as from any instrument that is networked (RS485). The following examples are for the HygroLab. The same commands can be used for other instruments with similar capabilities by substituting the appropriate values for the product ID and address.

1) Read %RH and temperature from the HygroLab 2 with address 01 (2 probes), directly connected to the PC COM port (RS232):

data request: {b01RDD}+CR

answer: {b01RDD 0025.01;0016.89;0024.57;0019.84;----.--;----.--;----.--;----.--;#C

note: 0025.01: humidity probe 1 (%RH)
0016.89: temperature probe 1 (°C as per instrument configuration)
0024.57: humidity probe 2 (%RH), etc.

2) read %RH, temperature and computed parameter

data request: {b01RDD0;}+CR

answer: {b01RDD 0025.90;0015.82;-003.69;0024.47;0019.88;-001.00;----.--;----.--;----.---S

note: 0025.90: humidity probe 1 (%RH)
0015.82: temperature probe 1 (°C as per instrument configuration)
-003.69: dew point probe 1 (°C as per instrument configuration), etc.

3) Read %RH and temperature from the HygroFlex 3 (1 probe) with address 00, connected by RS 485 to the HygroLab 2:

data request |{M00RDD} +CR note the | character (ASCII 124) preceding the command

answer: {M00RDD 0020.41;0019.87;----.--;----.--;#E

note: 0020.41: humidity probe 1 (%RH)
0019.87: temperature probe 1 (°C as per instrument configuration), etc.

Appendix 7: Water Activity Measurement

Water Activity: definition and applications

Definitions

The moisture content of a product can be defined as the percentage weight of water in relation to the dry weight of the product.

Products in which moisture can be present can be classified in two categories: hygroscopic and non hygroscopic. Examples of hygroscopic materials are salts, vegetal fibers, most metal oxides, many polymers, etc. Examples of non hygroscopic products are metal powders, glass granules, etc.

Regarding the moisture content of a product, we define static equilibrium as a set of conditions under which the product does not exchange any moisture with its environment. Under conditions of static equilibrium, the moisture content of a hygroscopic product depends on the nature of the product and also on the two following factors:

- (a) the partial pressure of water vapor in the immediate environment of the product
- (b) the temperature of the product

If the moisture content of a product is not dependent on both these factors, then the product is not hygroscopic.

Hygroscopic products may absorb water in different ways: sorption with formation of a hydrate, binding by surface energy, diffusion of water molecules in the material structure, capillary condensation, formation of a solution, etc. Depending on the absorption process, water is bound to the product with more or less strength. Moisture content can include both an immobilized part (e.g. water of hydration) and an active part.

Water activity A_w (or equilibrium relative humidity %ERH) measures the vapor pressure generated by the moisture present in a hygroscopic product.

$A_w = p / p_s$ and $\%ERH = 100 \times A_w$, where:

- p : partial pressure of water vapor at the surface of the product
- p_s : saturation pressure, or the partial pressure of water vapor above pure water at the product temperature

Water activity reflects the active part of moisture content or the part which, under normal circumstances, can be exchanged between the product and its environment.

Water activity is usually defined under static conditions of equilibrium. Under such conditions, the partial pressure of water vapor (p) at the surface of the product is equal to the partial pressure of water vapor in the immediate environment of the product. Any exchange of moisture between the product and its environment is driven by a difference between these two partial pressures.

Finally, water vapor can also be present in a gas or gas mixture. The relative humidity of a gas is defined as $\%RH = 100 \times p/p_s$, where (p) is the partial pressure of the water vapor present in the gas mixture and (p_s) is the saturation pressure, or the partial pressure of water vapor above pure water at the temperature of the gas.

Aw and Temperature

Both water activity (materials) and relative humidity (gases) are referred to the saturation pressure (ps) or partial pressure of water vapor above pure water:

$$A_w = p / p_s$$

$$\%RH = 100 \times p/p_s$$

The saturation pressure (p_s) is strongly dependent on temperature. At normal room temperature, (p_s) increases by about 6.2% for a 1°C increase in temperature. In an open environment that is not saturated with water vapor, the partial pressure of water vapor (p) does not change with temperature. In a closed environment, (p) changes proportionally to the °K temperature (°K temperature = °C temperature + 273.16). At normal room temperature, the change in (p) caused by a small change in the °C temperature is practically negligible. Because (p) does not change with temperature while (p_s) does, the relative humidity of a gas ($\%RH = 100 \times p/p_s$) is strongly temperature dependent. At 95 %RH and room temperature, an increase of 1°C results in a relative humidity decrease of about 6 %RH. At 50%RH, the same temperature increase causes relative humidity to decrease by about 3 %RH.

The water activity of most hygroscopic products is not as strongly dependent on temperature. At room conditions, research data typically shows that water activity varies only by roughly 0.0005 to 0.005 A_w (0.05 to 0.5 %RH) when temperature changes by 1°C.

This is explained by the fact that the partial pressure (p) at the surface of a hygroscopic product does vary with temperature. Above most hygroscopic products, the magnitude of the change in the partial pressure of water vapor (p) with temperature is similar (but not exactly equal) to the magnitude of the change of the saturation pressure (p_s) above pure water.

In summary, a change in temperature causes the partial pressure of water vapor above a hygroscopic product to change. At the same time, the partial pressure in the air above the product is practically unchanged. It follows that any change in the temperature of a hygroscopic product automatically causes the product to exchange moisture with the air (or gas) that surrounds it. Moisture is exchanged until the partial water vapor pressure at the surface of the product and in the air is equal. When measuring water activity, it is essential to keep temperature as constant as possible.

Applications

The active part of moisture content and, therefore, water activity, provide better information than the total moisture content regarding the micro-biological, chemical and enzymatic stability of perishable products such as foods and seeds. For similar reasons, water activity is equally relevant in the pharmaceutical industry where it provides useful information regarding the cohesion of tablets and pills, or the adherence of coatings. Water activity can be directly compared with the relative humidity of the ambient air to prevent dimensional changes in a product (paper, photographic film), to prevent hygroscopic powders (powdered sugar, salt) from caking or turning into a solid block, etc

Water activity can be used with some products (mostly synthetic products) as a means of indirectly measuring the total moisture content. This requires developing sorption isotherms to this purpose. Sorption isotherms are graphs that provide the relationship between water activity and moisture content at constant temperature. For most natural products, repeatable sorption isotherms cannot be reliably developed and water activity should be regarded as separate from moisture content.

Using the HygroLab 2 or 3 to Measure Water Activity

Both the HygroLab 2 and 3 can be used to measure the water activity of foods, pharmaceuticals and cosmetics. Both models can be configured to display humidity as water activity (1.00 aw = 100 %RH).

HygroLab 2:

When used as a stand alone unit, the HygroLab 2 can be used for conventional water activity measurements with up to 4 probes simultaneously.

Productivity in the laboratory can be greatly increased by using the HygroLab 2 together with a PC on which the HW3 or HW4 software has been installed. With this combination, the following functions are available (see note below):

- Accelerated water activity measurement (AwQuick mode): allows the measurement of most products in typically 5 minutes.
- Full equilibration measurement (AwE mode) with automatic detection of equilibrium conditions with optional recording of the entire equilibration process.
- Automatic capture to a text file of the measurements, including probe and product identification (both in the AwQuick and standard mode).

All HW3 water activity functions can be used simultaneously with up to 4 probes, with any or both measurement modes active.

Note: Access to the water activity functions of HW3 requires a product key. A product key can be purchased at the same time as the HygroLab or it can be purchased later. The product key works only with a specific instrument model and serial number. These functions are currently not available with HW4.

HygroLab 3:

Both the AwQuick and AwE measuring modes are already built into the HygroLab 3 and up to 4 probes can be used simultaneously.

When the HygroLab 3 is used together with a PC on which the HW3 software has been installed, the following additional functions are available:

- Recording of the entire equilibration process of a product sample (standard mode only).
- Automatic capture of measurements to text file (both in the AwQuick and standard mode).

No product key is required to access the water activity functions of HW3 with the HygroLab 3.

Note: the remaining pages of this addendum apply only to the HygroLab 3. The optional HW3 PC software features water activity functions that can be used with either the HygroLab 2 or the HygroLab 3. These functions are described in a separate manual provided on the HW3 CD ROM. These functions are currently not available with HW4.


Principle of Measurement

The standard (or conventional) method of measuring water activity consists in placing a sample of the product to be measured in a sealed container. The product sample slowly exchanges moisture with the air inside the sealed container until equilibrium is reached. Experience shows that even with a relatively small amount of product, it is the moisture contained in the product that ends up determining the humidity of the air, and not vice versa.

The equilibration process is monitored by measuring the humidity of the air above the product with a relative humidity sensor ($\%RH = 100 \times A_w$). Because temperature is an important factor when measuring water activity, the temperature of the air above the product, is also monitored.

By definition, water activity is equal to $\%RH / 100$ when equilibrium has been reached. At that time, the product no longer interchanges moisture with the surrounding air.

The following below uses the Rotronic AW-DIO and AwVC-DIO water activity probe as an example.

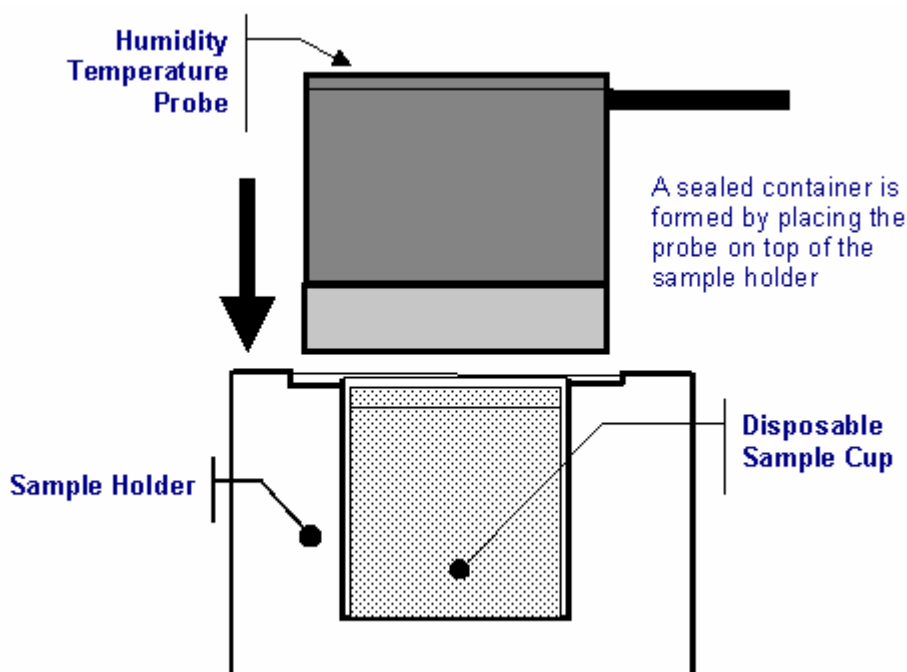

General Recommendations

1. Preparation of the product samples


Both the AW-DIO and the AwVC-DIO probes can be used with two different sample holders (WP-14: shallow / WP-40: deep), corresponding to two sizes of disposable sample cups (PS-14: shallow / PS-40: deep). Use the shallow sample cups (PS-14) and the shallow sample holder (WP-14) for calibrating the probe with the Rotronic certified humidity standards or with saturated salt solutions. In general, use the shallow cups for a liquid, a paste or a powder. Use the deep sample cups (PS-40) and the deep sample holder (WP-40) for products in bulk (large chunks).

The disposable sample cups serve two purposes:

- (a) provide a means storing product samples prior to measurements so that they can come to the same temperature as the probe
- (b) prevent contamination across samples.

Prior to measurements, fill a number of disposable sample cup with the products to be measured.

Usually, filling up to 1/3 of the cup is sufficient. To avoid soiling the probe, do not fill the cup above the fill line (see below).

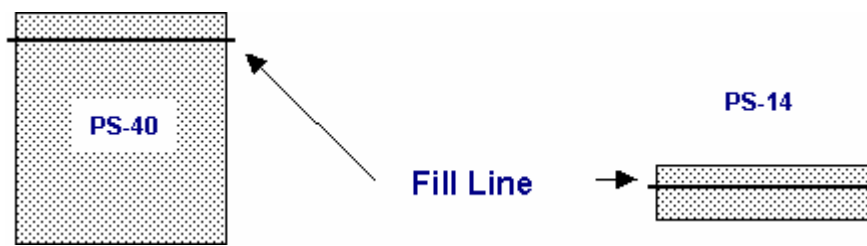

2. Other basic guidelines

To ensure accurate and repeatable results, please follow these basic guidelines:

- a) Measure water activity only in a temperature stable area. Do not measure on a bench that is located near a heater, an AC vent or an open window. Avoid direct exposure of the probe and/or product samples to sun light (heating effect).
- b) Prior to measurements, place each product sample in a disposable sample cup with the cover on. Place the cups in the same general area as the probe. Allow for sufficient time for the samples to come to the temperature of the probe (usually room temperature). A frequent mistake is to measure product samples which have been kept overnight in a refrigerator without first allowing them to come to room temperature. Another mistake is to measure hot samples coming straight from the manufacturing process. To measure water activity at a temperature other than room temperature, place both the probe and the product samples in an incubator set at the desired temperature.
- c) Avoid warming up the probe, the sample holder or the product sample by touching or holding for too long in your hand.
- d) When using either the AW-DIO or the AwVC-DIO probe, check that the probe is on (this is the default when powering up the HygroLab). When the probe is on, the red LED on the top of the probe should be flashing. The probe does not operate and the HygroLab does not get any signal when the LED is not flashing. If necessary, press once on the red button located on top of the probe.

Using the AwE Mode (HygroLab 3)

Note: be sure to configure the instrument to display humidity as Aw and temperature in °C.

In the AwE mode, the HygroLab 3 constantly monitors the humidity and temperature signals from the water activity probe. The HygroLab waits until the product sample has reached full equilibrium and automatically ends the measurement at that time.

The HygroLab 3 performs the following tasks:

- 1) the stability of the humidity signal is constantly evaluated
- 2) the stability of the temperature signal is constantly evaluated
- 3) the measurement is automatically ended as soon as both humidity and temperature reach equilibrium. At that time, the HygroLab makes a short beep and freezes the display. The RS 232 output keeps being updated.


The criteria used to determine that a product sample has reached equilibrium can be modified by the user (see Using the Standard Aw Mode).

The natural (or static) equilibration of most products typically requires from 45 to 60 minutes and can take as long as a couple of hours. Humidity changes at an extremely slow pace towards the end of the process. HW3 eliminates the tedious task of monitoring the measurements until equilibrium has been truly reached.

When ready to measure:


1. Press on the **ENTER** key on the keypad of the HygroLab 3. The display shows:

Note: the numerical values are for illustrative purposes only.

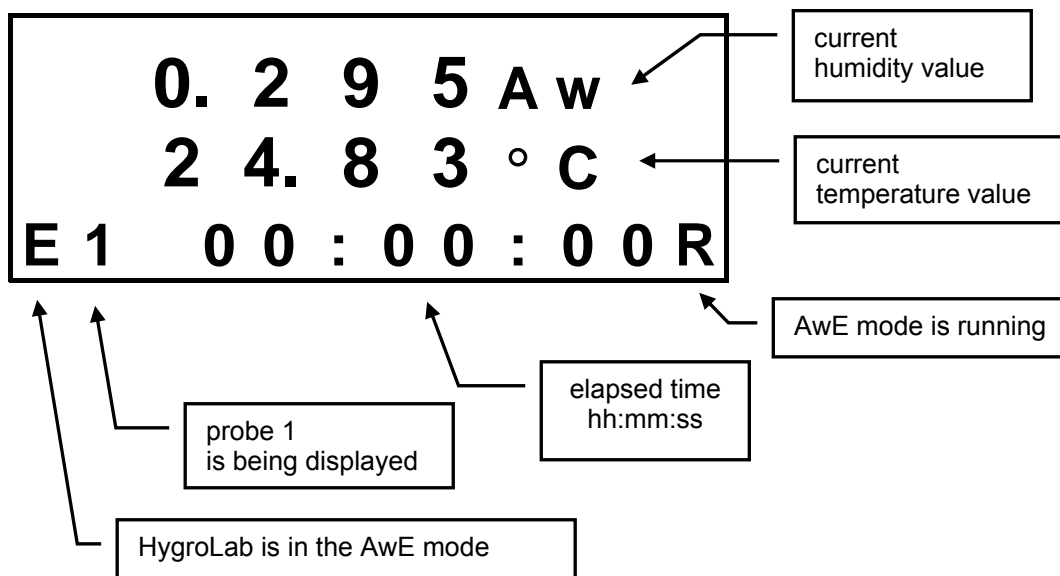


After 60 seconds, the trend indicators appear both for humidity and temperature. Two arrows mean stable signal, as per the settings defined when entering the AwE mode :



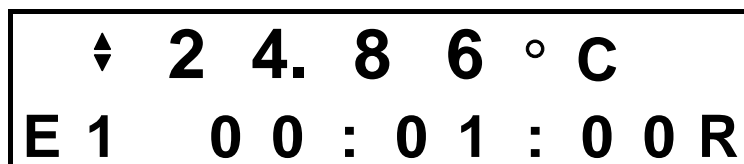

Upon reaching and holding equilibrium for a few minutes, the HygroLab beeps once and freezes the display. Black rectangles to the left of the display indicate that the measurement is ended.

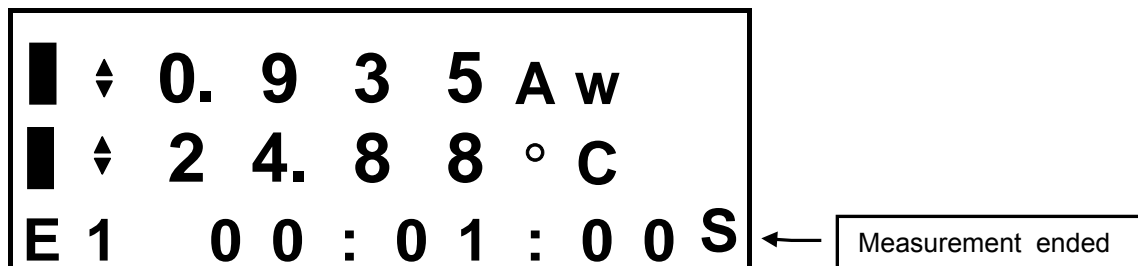

The AwE mode runs simultaneously for all probes connected to the HygroLab. When a probe is done, the HygroLab beeps once and freezes the display for that probe. The display can be manually switched to another probe from the keypad with the ▼ or ▲ arrow.

When the last probe is done, press the **ENTER** key to exit this screen and return to the standard mode.

Press **ENTER** again to start a new measurement in the AwE mode with all probes.

Important: unlike the equivalent function in the HW3 software, the HygroLab 3 does not allow to start and stop each probe individually. To avoid possible confusion, you may want to either disconnect or turn off any unused probe.

Using the AwQuick Mode (HygroLab 3)

In the AwQuick mode, the HygroLab 3 uses an algorithm to project the full equilibrium value (water activity) of the product sample. The measurement is automatically ended and typically requires about 5 minutes.

The HygroLab 3 performs the following tasks:

- 1) the value of the humidity signal is constantly monitored
- 2) the stability of the temperature signal is constantly monitored
- 3) After an initial period of time (dwell time), HW3 uses the humidity data to project the end value of the equilibration process (water activity). The measurement ends automatically as soon as the projected Aw value is stable. At that time, the HygroLab makes a short beep and freezes the display. The RS 232 output keeps being updated.

With the default dwell time of 4 minutes, measurements typically require about 5 minutes. The difference between the AwQuick mode and the Standard mode is typically 0.005 aw or less.


The value of the dwell time can be set by the user (see Mode). This value represents a trade off between speed of measurement and accuracy. Generally, a longer dwell time produces more accurate results but causes measurements to take longer.

The value of temperature is the average temperature during the measurement. The HygroLab displays a trend indicator next to the value of temperature. This is used to verify that temperature is stable during the measurement.


When ready to measure:

1. Press on the **ENTER** key on the keypad of the HygroLab 3. The display shows:


Note: numerical values are for illustration purposes only.

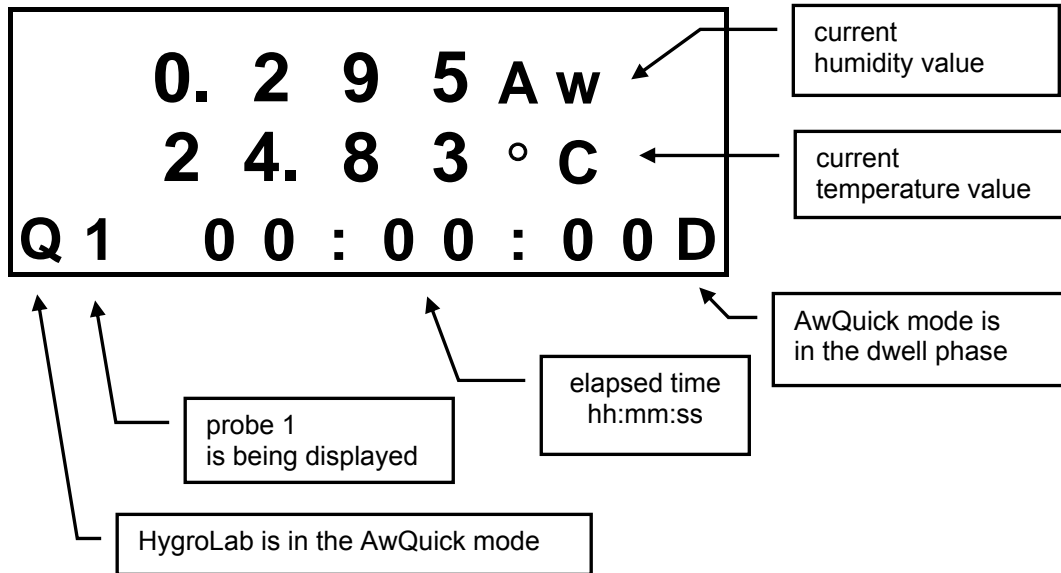


After 60 seconds, the trend indicator appears for temperature. Two arrows mean stable signal, as per the setting defined when entering the AwQuick mode:

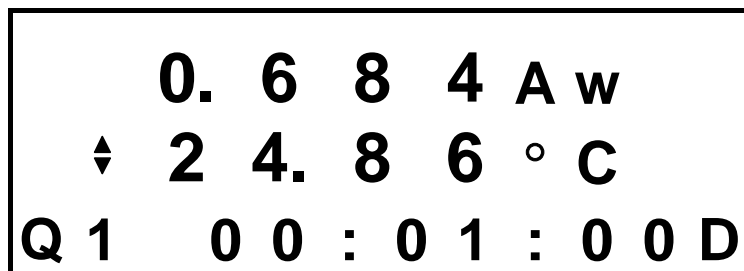

At the end of the dwell time, the display changes as follows:

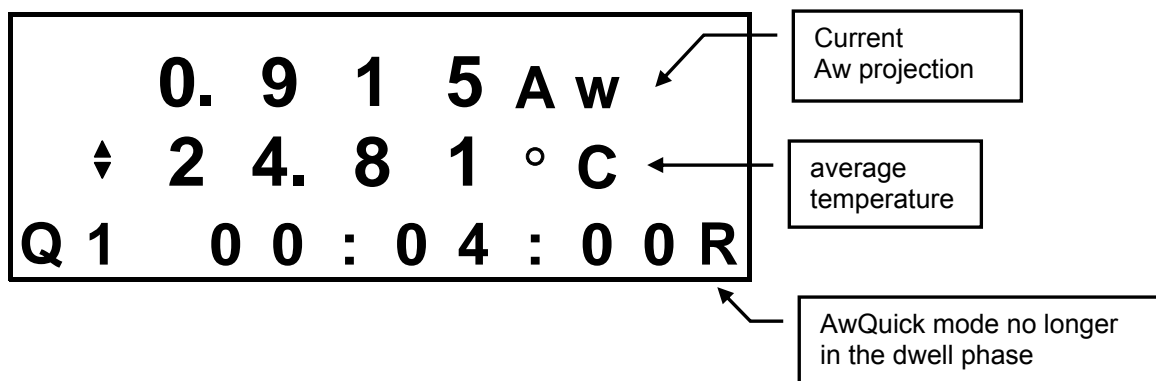

When the projected value of Aw is stable, the HygroLab beeps once and freezes the display. Black rectangles to the left of the display indicate that the measurement is ended.

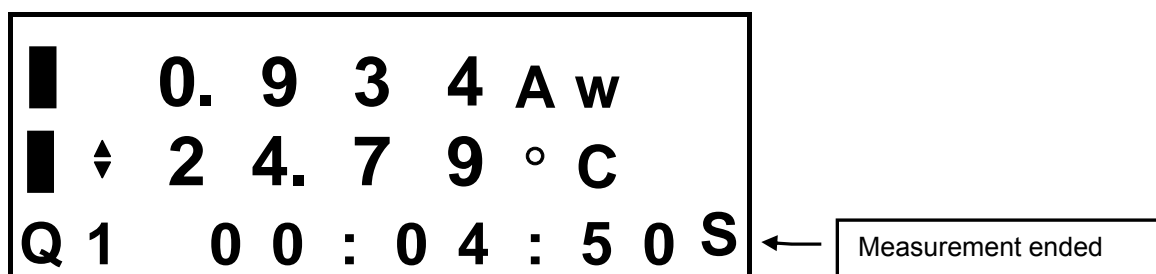

The AwQuick mode runs simultaneously for all probes connected to the HygroLab. When a probe is done, the HygroLab beeps once and freezes the display for that probe. The display can be manually switched to another probe from the keypad with the ▼ or ▲ arrow.

When the last probe is done, press the **ENTER** key to exit this screen and return to the standard mode.

Press **ENTER** again to start a new measurement in the AwQuick mode with all probes.

Important: unlike the equivalent function in the HW3 software, the HygroLab 3 does not allow to start and stop each probe individually. To avoid possible confusion, you may want to turn off any unused probe either by disconnecting or turning off the probe.

Appendix 8: Accessories for the HygroLab

Order Code	Description
HW3	HW3 software (CD ROM)
HW4-E	HW4 software, Standard Edition (single user)
HW4-P	HW4 Professional Edition (FDA/GAMP validation functions, multi user)
AC1623	Cable service connector (MTA) HygroLab 1 to PC
MOK-01-B5	Adapter cable for HygroClip probes with DAT05 connector. Cable length 1 meter (3.2 ft)
T7-03-B5	Adapter cable for HygroClip probes with T7 connector. Cable length 3 meter (9.8 ft)
MOK-03-WIN	Calibration cable for HygroClip probes with DAT05 connector. Terminated with a 25-pin SUB-D connector. Converter 25-pin to 9-pin is supplied. Cable length 3 meter (9.8 ft). 9VDC adapter may be required (consult your local distributor)
T7-03-WIN	Calibration cable for HygroClip probes with T7 connector. Terminated with a 25-pin SUB-D connector. Converter 25-pin to 9-pin is supplied. Cable length 3 meter (9.8 ft). 9VDC adapter may be required (consult your local distributor)
EA00-SCS	0%RH humidity std, SCS cert., pack of 5
EA05-SCS	5%RH humidity std, SCS cert., pack of 5
EA10-SCS	10%RH humidity std, SCS cert., pack of 5
EA11-SCS	11%RH humidity std, SCS cert., pack of 5
EA20-SCS	20%RH humidity std, SCS cert., pack of 5
EA35-SCS	35%RH humidity std, SCS cert., pack of 5
EA50-SCS	50%RH humidity std, SCS cert., pack of 5
EA65-SCS	65%RH humidity std, SCS cert., pack of 5
EA75-SCS	75%RH humidity std, SCS cert., pack of 5
EA80-SCS	80%RH humidity std, SCS cert., pack of 5
EA95-SCS	95%RH humidity std, SCS cert., pack of 5
ER-15	calibration device for 15mm dia. Probes
ER-05	calibration device for 5mm dia. Probes
EM-15	calibration device for type 'IE' probes